


**XIEM  
EX**

# Extraction Balloon Catheter

<Short Wire Type>

Made in JAPAN

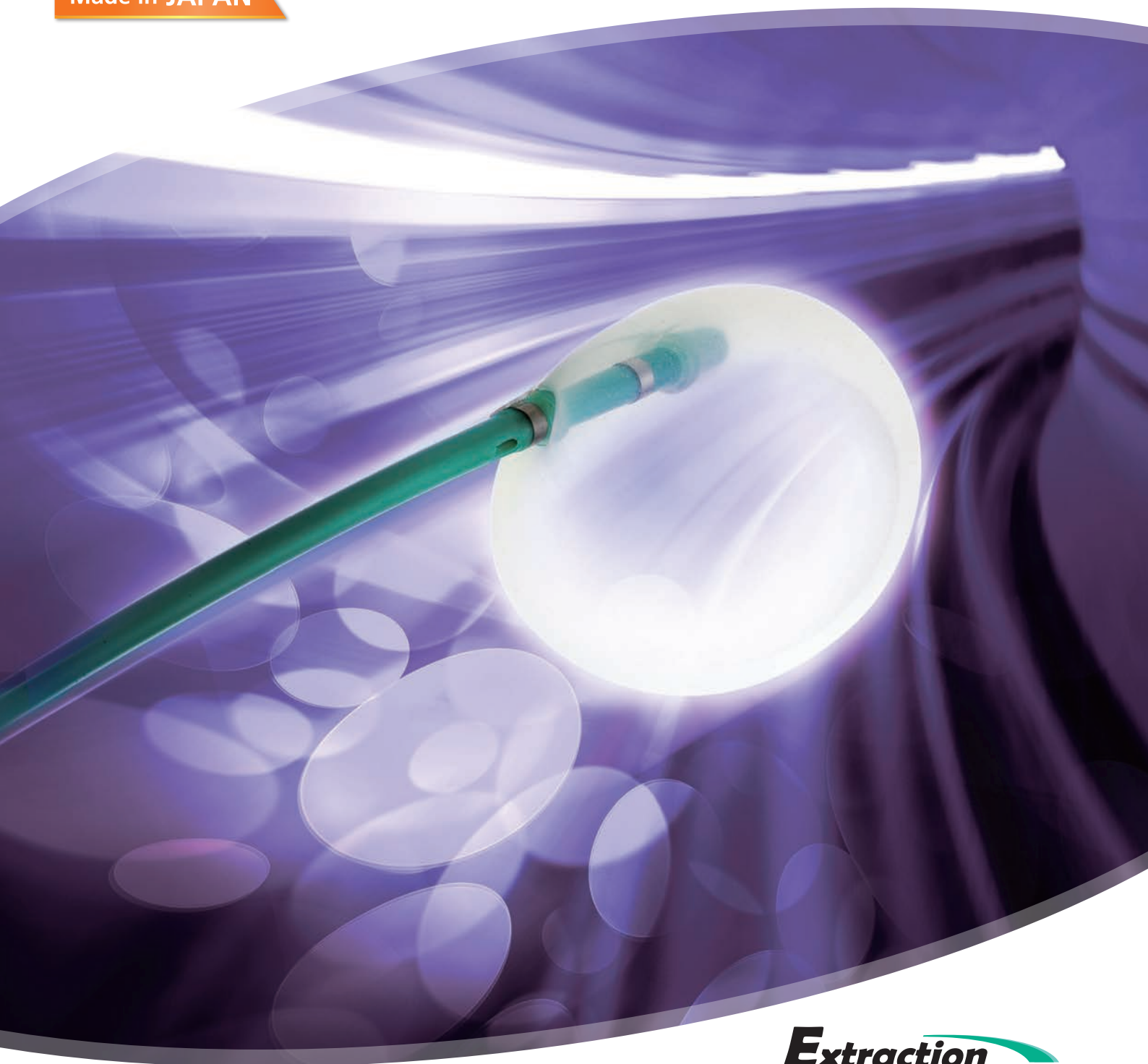

**Extraction**  
**Balloon Catheter**

ZEON MEDICAL INC.

CE 0197


Dedicated catheter that combines flexibility and pushability and aims to approach stones with smooth and simple repetition.



## Experience the performance of a unique balloon shape (P type balloon)

The P type balloon prevents the deformation of its shape even when it is subjected to external pressure from bile stones and the bile duct.

For higher performance removing biliary stones, sludge and sand, the P type balloon stays in close contact with the bile duct wall and it can enter dead spaces in the lower part of the bile duct.

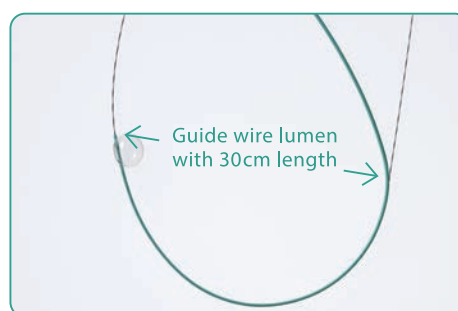

Product Code	Catheter size		Balloon shape & size			Compatible channel inner diameter (mm)	Injection hole	Number per package
	Outer diameter (F / mm)	Effective length (cm)	Shape	Expansion diameter (mm)	Volume (mL)			
EXP71820PM(E)	7-5.5 / 2.4-1.85	200	P type (Offset)	11 - 15 - 18	5	≥ 3.7	Proximal	1 pce
EXP71820DM(E)							Distal	1 pce

Applicable guide wire diameter is 0.88 mm (0.035 inch)

Guide wire is compatible with short wire lumen and this product cannot be directly inserted or inserted over the wire.

**Note**

- Ensure that this product is used by medical doctors or healthcare professionals who are taken suitable training for clinical endoscopic procedures.
- Be sure to always read the attached IFU when using this product.
- The specification and appearance of this product may change without prior notice.

Manufactured by:

**ZEON MEDICAL INC.**

1-6-2 Marunouchi, Chiyoda-ku, Tokyo 100-0005, Japan  
Phone: +81-3-3216-1266 FAX: +81-3-3216-1270

XEMEX is the brand by ZEON MEDICAL INC.  
Copyright ©2019 ZEON MEDICAL INC. All rights reserved.  
Ver. Jan. 2019 (E)EX-014 0119005(OC01)

Agent: