

ES DILATOR -Soft type-

Universal catheter dilator



ES Dilator
Effective Dilator for Interventional EUS

ZEON MEDICAL INC.

SHARPNESS

The tip was made thinner to reduce the height difference according to the size of the guidewire.




HIGH PUSHABILITY

Fluoride resin, which combines the flexibility and pushability of the tube material, enables the up-table operation, smooth surface insertion, and smooth insertion.

Good visibility

An Imaging ring located 10 mm from the tip can be used to identify the tip position under fluoroscopy. The Imaging ring is inside the catheter, so the outside of the catheter is smooth.

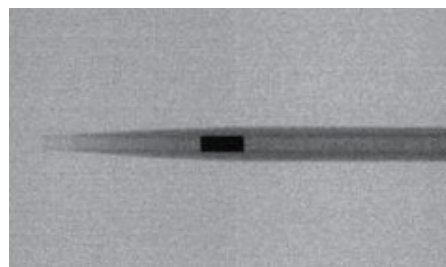

SPEC

P / N	Catheter outer diameter	Catheter effective length	Applicable guide wire
DC7P180E(E)	2.45mm(7F)	1800mm	0.46mm (0.018 inch)
DC7F180E(E)	2.45mm(7F)	1800mm	0.64mm(0.025 inch)
DC7R180E(E)	2.45mm(7F)	1800mm	0.89mm(0.035 inch)

attention

• Before using product, be sure to read the attached IFU.

manufacturer

ZEON MEDICAL INC.

[H O] Shin-Marunouchi Center Building, 1-6-2 Marunouchi, Chiyoda-ku, Tokyo 100-0005, Japan
Phone: +81-3-3216-1265
URL : <https://www.zeonmedical.co.jp/e/>

XEMEX is a house mark of Zeon Medical Inc. and is registered in Japan and other countries.